PHYSIOTHERAPY & REHAB


1. Interferential, TENS, Electro-stimulator, Combination Therapy

PhySys= Ultimate Combination Therapy with 12" colour touch screen from Zimmer Germany



All low & medium frequency currents & multi frequency U/s

2 Independent Channels- Parallel - Separate- Alternate with Channel Inversion

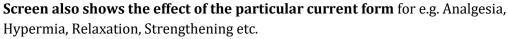



Current for Sports Rehab- General Muscle Strengthening, Improvement of muscle endurance, Strengthening of fast-twitch fibre muscles, endurance training, fast-twitch fibre muscle training, Explosive strength training, Sequence training- warm up- train-relaxation etc.

Special Treatment- Vein Strengthening, Psoriasis, Electro-Sleep, Lymph Oedema and Hyperhidrosis etc

Currents for Urology/ Gynaecology- Pelvic floor stimulation, Incontinence Treatment


Paralysis diagnosis & Paralysis treatment

Lange Test, Neofaradic test, Fischgold test

Rheobase, Chronaxie, accommodation quotient and I/T curve plotting

Ultrasound:

- Sono Swing Choice of penetration depth (Zimmer Patent)
- ➤ Choice of Mechanical & Thermal effect

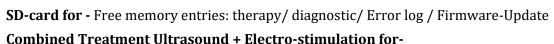

Hyperaemia, metabolism stimulation, pain control, muscle relaxation, trigger point treatment



SOLEO LINE- COMBO, Electrotherapy &Ultrasound 8" colour touch screen from Zimmer



Soleo SonoStim-Combination Therapy = Soleo Galva + Soleo Sono

Dual Channel for All medium & Low frequency currents & multi frequency Ultrasound All low & medium frequency currents & multi frequency Ultrasound Therapy



Screen also shows the effect of the particular current form for e.g. Analgesia, Hypermia, Relaxation, Strengthening etc.



Current for Sports Rehab- General Muscle Strengthening, Improvement of muscle endurance, Strengthening of fast-twitch fibre muscles, endurance training, fast-twitch fibre muscle training, Explosive strength training, Sequence training- warm up- train-



2. Ultrasound Therapy

SOLEO SONO Ultrasound 8" colour touch screen from Zimmer



Soleo Sono Ultrasound Therapy unit - Multi Frequency U/s Head offering choice of penetration



Ultrasound: Zimmer Patent 5 Choice of penetration depth



Ultrasound Frequency of 0.8 MHz for deep and 2.4 MHz for superficial effects and possibility of combining two frequencies in single treatment for multiple choice of penetration depth.

Ultrasound head has audio & Visual contact feedback with coupling % displayed on screen as well

Ultrasound head offers under water treatment option on screen SD Card for memory, software upgrade and service facility



Sono One Ultrasonic therapy unit from Zimmer Medizin Germany Features

- •Multi frequency, 1 MHz: deep effect / 3 MHz: superficial effect
- •Transducer head with visual contact indicator.
- •Colour touch screen monitor for better control.
- •Continuous & Pulsed: Pulse frequencies 10, 20, 50, & 100 Hz

3. Multiwavelength High Power Laser Therapy Unit



OPTON PRO 25W, 15W & 10W, therapeutic Laser for deep penetration & instant effect from Zimmer Germany

12" Colour touch screen

Operating mode - Continuous, Serial pulse, Burst, Single pulse

Integrated Visual analogue scale- VAS provides information about the patient's pain perception.

Thermal threshold test assists with calculation of the required dosage in combination with the skin type and the required penetration depth & prevents overdosing



3 Wave Lengths 810nm, 980nm, 1064nm

Therapy with 3 wavelengths enables stimulation of different target sites: Pain receptors at the surface as well as deeper tissue structures are stimulated at the same time. Wavelength of 1064 nm having lowest scatter in the tissue promotes penetration of the radiation into the deep layers of tissue.





2 Wave Lengths 810nm & 980nm

Therapy Assistant Mode-



It's Easy – Select body region > Select skin type > Determine therapy area By determining the therapy area, the course of the disease and the skin type, the following information is automatically saved:

- ✓ Program parameters
- ✓ Amount of energy to be applied
- ✓ Therapy time



Laser applicator with manual switch - the laser light is now emitted continuously, without any further action needed

Integrated Skin temperature measurement sensor enables the skin temperature to be monitored during the therapy.

4. Cryoflow Cryotherapy Unit



Cryo 7 & Cryo 6 – The new standard in cold air therapy, **from Zimmer Germany**, **"World largest Manufacturer of Cryo"**

Ÿ Brilliant, high-resolution, 10-inch touch display.

Ÿ Individual programmes can be stored.

Ÿ Optional supporting arm facilitates hands-free operation.

Ÿ Glass shelf plate for use as a tray for additional devices.

-35 °C Cold Air Cryotherapy for Pain and Inflammatory Treatment

Indications:

Painful acute and chronic conditions of the musculoskeletal system

Pain relief in rheumatic disorders

Neurological disorders, reduction in spasticity

Prevention of edema and hematomas

Targeted treatment of muscle trigger points in combination with stretching. Vascular training,

As pre-treatment of CPD (complex physical decongestion) in chronic lymph edema **Glass touch pad. No Consumables**


5. Focussed Shockwave Therapy Unit



enShock -Piezoelectrically <mark>Focused shock wave</mark> Therapy

(ESWT) with focus on healing from Zimmer Germany

- > Long-term pain relief
- Positive effect on inflammatory reactions in musculoskeletal discomfort/conditions,
- Excellent effects on bone healing, especially in non-healing fractures and pseudo-arthrosis
- **Elimination of calcification in tendons**

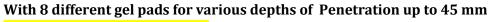

11" Colour touch screen

6. Radial Shockwave Therapy Unit



12" Colour touch screen EnPuls Pro from Zimmer Germany

Radial Shock Wave Therapy unit with **electro-magnetic** shock generation & Colour touch screen for ease of operation

Special Features

- **✓** Ballistic radial shockwave Electromagnetic generator
- ✓ Efficient down to 35 mm depth
- ✓ Radial tissue dispersion
- √ High mass projectile
- · Ports for two hand pieces
- · Energy levels from 60 to 185 mJ (equivalent 1-5 bar)
- · Pulse frequencies up to 22 Hz & Burst mode for Trigger Points

Less pain, same energy – enPuls soft shot technology Medical effects

- ✓ Immediate, sustaining pain relief
- ✓ Increased range of mobility
- ✓ Controlled normalized muscle tone
- √ Change of membrane permeability



8" Colour touch screen enPuls 2 from M/S Zimmer Germany

Portable Radial Shock Wave Therapy unit with **electro-magnetic** shock generation & Colour touch screen for ease of operation

Less pain, same energy - enPuls soft shot technology

- > Energy levels from 60 to 185 mJ (equivalent 1-5 bar)
- > Pulse frequencies from 1 to 22 Hz
- > Indication menu with preset treatment protocols

7. High Energy Inductive Therapy-Magneto-therapy @ 3 Tesla Power

Power of MRI for Physiotherapists EmField Therapy Unit @ 3 Tesla from Zimmer Germany



Exposure to natural energy fields is a requirement of human life. These energy fields interact with our cells, and are necessary for homeostasis & maintaining our bodies at an optimal level.

Zimmer Medizin Systeme GmbH Germany introduces High Energy Inductive Therapy unit "emFieldPro" for healing in a natural way by Earth's magnetic field (**EmField Therapy**)

The **emFieldPro** generates a magnetic field of **3 Tesla**, which is about 600 times stronger than a normal magnet bar.

This strong magnetic field stimulates nerve cells, muscles, and blood vessels consequently as this application of energy results in a physiological change and/or stimulation which is used to achieve therapeutic effects.



EmField Therapy ADVANTAGES:

- Non invasive.
- Totally safe procedure applied from distance without contact
- Easy and fast application even over clothing
- Short treatment duration
- Low operation cost without any consumables



EmField Therapy offers a safer and more comfortable alternative for urinary incontinence to prior treatments.


EmField Therapy improves sports performance, and simply helps to maintain good health.

EmField Therapy stimulates muscles, connective tissues, intestines, tendons and cartilage, the brain and peripheral nerve sites.



INDICATIONS:

Spinal Injuries/Disorders, Chronic/Acute Lumbago Chronic/Acute Ischialgia, Spondylitis, Nerve Disorders, Damaged Peripheral Nerves, Degenerative Arthritis, Rheumatoid Arthritis Cervical Syndrome, Frozen Shoulder, Prostate Pain

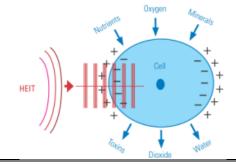

8. Hands Free TECAR Therapy



ThermoTK, operating at frequencies of **460 kHz and 540 kHz,** lies in the lower region of the high-frequency treatment.

ThermoTK is a long-wave therapy device.

Indications: • Pain Management • Acute –and Chronic Pain • Inflammations • Improves blood circulation • Muscle relaxation • Hypertonic muscle relaxation • Post-training regeneration • Trigger points

Treatment Modes:

Active: (optional with 2 applicators)

- Continuous mode with 460 or 540 kHz
- Pulsed mode with 460 or 540 kHz
- Alternation between 460 and 540 kHz in continuous mode

Passive, hands free (always resistive mode)

9. Inductive Shortwave Diathermy



Thermo Pro – Modern Shortwave Therapy Unit - From Zimmer Germany Solid State Shortwave Diathermy for trouble free operation.
On screen coupling display ensures optimum energy transfer Pre programmed indication menu as per body area and facility of storing own programs in memory.

Inductive Coil based Circuplode for radiations free Treatment. Medical effectiveness

- Hyperaemia, Muscle relaxation, Metabolism boosting, Re-absorption promotion

10. Cervical & Lumbar Traction System

LUMBER & CERVICAL TRACTION



clTrac - Cervical & Lumbar Traction System A safe alternative to pain medication:

Traction therapy with clTrac is highly effective for intervertebral disc problems and discomfort in the back and neck muscles: The new digital traction device from Zimmer expands the joints of the cervical spine, lumbar spine, and the extremities (lumbar) with metered tensile force.

This provides immediate relief to irritated nerves and the surrounding tissue can regenerate and heal.

Indication-based with software-controlled setting parameters

You can regulate all parameters in small increments,

1 to 99 minutes and 0.1 Kg to 90Kg and switch between digital and mechanical control.

Touch screen display for ease of operation.

11. Decompression Therapy System

DECOMPRESSION THERAPY

GRAVITY ASSISTED SUPINE

ANTALGIC POSITIONING

EXTENSION TREATMENT

FLEXION TREATMENT

Lumbar...Prone...Supine...Flexion...Inversion...Vibration! Side Lying, Cervical

Gravity Assisted Supine - ONLY Decompression Table featuring a bolster-less leg elevation with pelvic tilt which pre-tensions posterior elements.

Antalgic Positioning: Allow for the safe and comfortable treatment of acute and antalgic patients who exhibit non-weight bearing relief postures.

Extension and flexion in the lumbar spine

Axial lateral pulls, lumbar & Cervical for lateral discs Bolster less leg elevation

Treat Thoracic problems... unload the "entire spine" equally by un-weighting with inversion.

With a centre-section pelvic drop we can now address SI joint fixation/asymmetries as well as the common PI/AS

Powered Caudal Sections: Controlled by 3-position hand switch (elevation, thoracic elevation and lumbar elevation/flexion

The Passive Tensioning Reactive Orthopedic Spring controls X-Axis Table motion. Enhance decompression without added force. Tension can be matched to patient morphology.

Decompression to the Cervical Spine- mid-section elevation allows full support to the lumbar spine while lowering Thoracic kyphosis.

Enables a full range of cervical angulations

Allow cervical discs to migrate centrally without an increase in muscular tension.

CARPAL, HIP & WRIST Traction KIT